

MEDI CAPTURE

™ Simple  
Connected  
Powerful



NEW  
RELEASE

## MVR Lite

The MVR Lite is the simplest in the series of the MVR medical video recorders from MediCapture. It is developed specifically for outpatient surgical facilities that are looking for a video recorder that is easy to use and low cost, but still comes with many of the advanced features of a higher-end recorder that are critical when recording surgeries today.



MVR Remote App<sup>3</sup>  
(with Connect Upgrade Package<sup>1</sup>)



### Compact and Powerful

Compact, rugged, and so powerful, the MVR Lite records in full HD at 60 fps and offers unique features, such as dual input for two-camera operation.



### Smart Workflow

MVR Lite offers a unique open user interface for easy control, including touchscreen operation and instant access to archives. Managing recorded data has never been so easy!



### Upgradable for Evolving Needs

MVR Lite is built to “evolve” with the needs of medical facilities with various upgrade options, including 4K recording or the Connect Upgrade Package that offers network storage and tablet control.<sup>3</sup>

Medical Imaging Made Easy™

# Features



## Multi-Input Recording

Directly toggle between two inputs and perform picture-in-picture recordings.



## MVR Manager

Browse and review all patient files, annotate, zoom, report.



## MVR Secure

Protect your files from hacking or power outages. Includes user authentication and historical logs for data protection.




## Touchscreen Operation<sup>4</sup>

Connect the MVR Lite to a USB touchscreen monitor for even easier operation.

# Upgrade Packages

The MVR Lite Upgrade Packages can be installed at any time by simply using an assigned activation key provided by MediCapture.

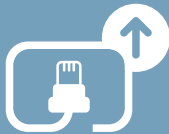

## Connect<sup>1</sup>

Provides network storage, tablet control through the MVR Remote App<sup>3</sup>, universal printing and Administration access.



## 4K Recording<sup>1,2</sup>

Records delicate surgeries in stunning 4K (2160p30/25).

<sup>1</sup> Option

<sup>2</sup> Check camera compatibility with MediCapture prior activation

<sup>3</sup> Compatible with Android tablets and mobile phones

<sup>4</sup> List of compatible Touchscreens on our website

<sup>5</sup> Subject to change without notice

# MVR Lite Specs<sup>5</sup>

<b>Device type</b>	<b>HD Medical Video Recorder</b> Records high-definition video and still images from medical video sources to external USB drive and network server (with Connect Upgrade Package <sup>1</sup> )
<b>Case/Buttons</b>	Stainless steel, fluid resistant
<b>User interface</b>	Overlay on external monitor via Touchscreen / Mouse / Touch-Pad.
<b>Video Inputs</b>	HDMI, DVI (with adapter)
<b>Video Outputs</b>	HDMI, DVI (with adapter)
<b>Video Formats</b>	MPEG4. H.264 / H.265
<b>Input video resolutions</b>	1920x1200p60/50, 1080p60/50/30/25, 1080i60/50, 720p60/50, 1280x1024p60/50 Input 1 for 4K also: 2160p30/25 <sup>1</sup>
<b>Image Format</b>	JPG / PNG
<b>Sound Recording</b>	Mic by Stereo-Headset (4 poles) / USB microphone
<b>External Storage</b>	USB flash drive, USB hard drive (FAT32/NTFS/exFAT)
<b>USB Support</b>	USB3.0, SuperSpeed, 1 front connector, 2 rear connectors
<b>Footswitch/Trigger</b>	3.5mm mini jack. Compatible with standard footswitches (not included) and closed-contacts camera triggers
<b>Remote Control</b>	IR remote control (included), MVR Remote App <sup>3</sup> (with Connect Upgrade Package <sup>1</sup> )
<b>LAN<sup>1</sup></b>	USB-LAN Adapter required (NA400USB) RJ-45, Ethernet 10/100/1000BASE-T
<b>Mechanical</b>	Size: 160 x 140 x 40mm / 6.3 x 5.5 x 1.6" Weight: 0.87 kg / 1lb 14,7oz
<b>Environmental</b>	Storage: -40° to +85° C Operating: -20° to +40°C
<b>Power</b>	AC 100-240V, 50/60Hz, 0.4-0.2A, 3-prong AC jack
<b>Certifications</b>	CE 93/42/EEC ANSI/AAMI ES60601-1 (2005(R)2012 + A:2012, C1:2009(R)2012 + A2:2010(R)2012) (EN 60601-1 ed 3.1), CAN/CSA C22.2 No.60601-1-1:2014-Edition 3, FCC Part 15, Class 1 Equipment according to US-FDA, RoHS3 / ISO 13485: 2016
<b>Order number</b>	MVR400 MVR Lite recorder, USB flash drive, HDMI cable, HDMI to DVI adapter, split remote cable adapter, power cable, IR remote control, User Guide (English)
<b>Order information for optional upgrade packages:</b>	AK4K: 4K Recording Upgrade Package NA400USB + AKNETWORK: Connect Upgrade Package

MVR Remote App<sup>3</sup>  
(with Connect Upgrade Package<sup>1</sup>)



## MediCapture, Inc. – North America

2250 Hickory Road, Suite 200  
Plymouth Meeting, PA 19462, USA

+1-888-922-7887

www.medicapture.com  
sales@medicapture.com

## MediCapture Europe B.V.

Kraanpoort 10  
6041 EG Roermond  
The Netherlands

+31 475 214 305

www.medicapture.com  
sales.europe@medicapture.com

## MediCapture APAC Inc.

3FL-6, No. 432, Section 1, Keelung Road,  
Xinyi District 110, Taipei City,  
Taiwan, Republic Of China

+886- 2- 2345- 3668

www.medicapture.com  
sales.asiapacific@medicapture.com